

exofin[®] PRECISION PEN

High Viscosity Topical Skin Adhesive



Design Innovation for Precise Application

Every patient is different. Every procedure is unique. For effective wound closure and optimal healing, the skin adhesive needs to be applied with precision — regardless of the wound size or location. To achieve the level of control necessary to close even the most delicate procedures, the adhesive applicator should be natural to hold, safe to operate, and easy to use.



With its glass-free, ergonomic design, multiple application options, variable flow rate, and industry leading highest viscosity, the exofin® Precision Pen enables medical professionals to apply the industry's most advanced skin closure adhesive safely and efficiently — with adaptable precision.


Along with exceptional comfort and control, the Precision Pen also features the fastest drying adhesive and a 14-day microbial barrier with long-lasting protection.

exofin[®] PRECISION PEN

THE PERFECT LINE EVERY TIME

FEATURES & BENEFITS

- Ergonomic design with grooved, non-slip surface
- Pressure-controlled, variable adhesive flow rate
- See-through activation chamber for easy viewing
- 152% more usable adhesive than Dermabond Advanced[®]
- Highest viscosity
- Fastest drying adhesive
- 14-day microbial barrier

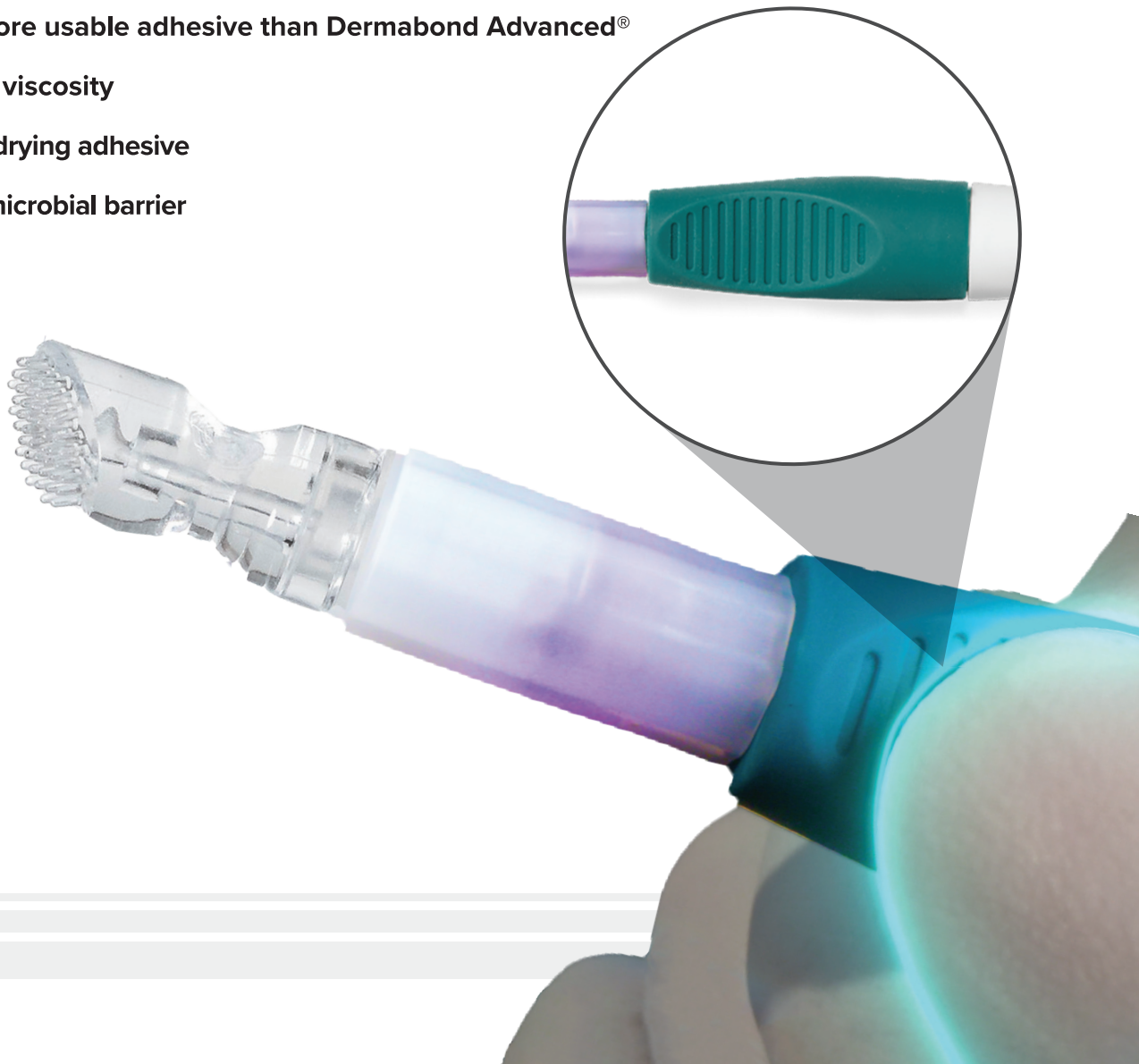


Performance Comparison

TEST METHOD	exofin [®] PRECISION PEN	DERMABOND [®] ADVANCED [®]	LIQUIDBAND [®] EXCEED [™]	LIQUIDBAND [®] RAPID [™]	LIQUIDBAND [®] FLOW CONTROL/ SWIFTSET [™]
510K Number	K212246	K152096	K151182	K183570	K122446
Compositional Analysis	100% 2-OCA	100% -OCA	100% 2-OCA	100% 2-OCA	100% n-BCA
Tensile Strength (MPa)	0.227	0.224	0.167	0.147	0.188
Wound Closure Strength (N)	17.24	15.97	8.82	13.46	10.59
Heat of Polymerization ($\Delta^{\circ}\text{C}$)	5.49	9.02	2.33	7.5	3.17
Cure Time (sec)	56	91	*	47	*
Viscosity	351 cP	262 cP	63 cP	40 cP	5 cP
Adhesive Expression (g)	0.634	0.252	0.2088	0.2999	0.2117
Flexibility [*] (# of cycles to failure)	300	117	277	145	1


* application dependent




THREE APPLICATIONS FOR VERSATILE CLOSURE

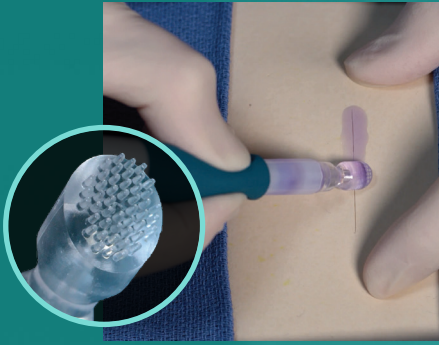

**Micro-Bristle Brush for
Wide Coverage**

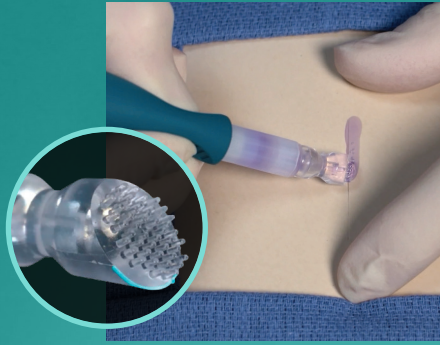

**Inverted Brush for
Narrow Lines**

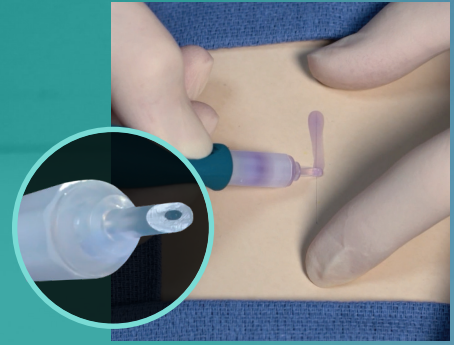

**Angled Precision Tip
for Micro Lines**

The exofin Precision Pen is the only pen-style applicator on the market with a variable flow rate and an adjustable tip that allows the user to control the adhesive width. Every aspect of this innovative applicator is designed to improve clinician control and patient outcomes.

To Place an Order: Email: customerservice@chemencemedical.com

Contact Info: Phone: 844-MED-GLUE (633-4583)



All Contents © Copyright 2021 Chemence® Medical, Inc. All rights reserved.

No modification or use, or further reproduction or copying, of any of the content is authorized without the express consent of Chemence® Medical, Inc.